

# PAL®

## Power-Assisted Liposuction


The PAL® System from MicroAire Aesthetics is a patented surgical instrument system used for fat harvesting and body contouring.

**PAL® has very real benefits over other liposuction products:**

1. It uses no heat or ultrasound energy, so the fat cells are just as healthy and viable as cells harvested using manual liposuction.
2. It uses reciprocating motion to pass easily through tissue, especially fibrous areas, which makes procedures faster, and results in less bruising and bleeding.

**The PAL® System consists of four components:**

- The **PAL-650 handpiece** has a thumb-throttle to control the rate of reciprocation. It is engineered for ultra-low vibration, comfort, efficiency and precision control.
- The **NEW 5020 Electric Console** provides reliable power to the handpiece, utilizing a touch-screen interface that can be used to set the instrument's max speed. The console accepts two instruments simultaneously, including any of MicroAire's small-bone orthopedic power instruments. The low-flow irrigation pump is intended for orthopedic use, not for tumescent infiltration. The 5020 Console also meets CE requirements for reduced electrical interference that went into effect in 2013. The 5020 does accept older PAL® handpieces, but a new 5006-PAL Smart Cable is required.
- The **5006-PAL Smart Cable** connects the handpiece to the console and tells the console which type of instrument is connected. Cable is 3.66 meters long.
- **MicroAire Aesthetics Cannulas** are offered in a wide variety of sizes and port configurations. All cannulas are made of high-quality, electro-polished stainless steel.

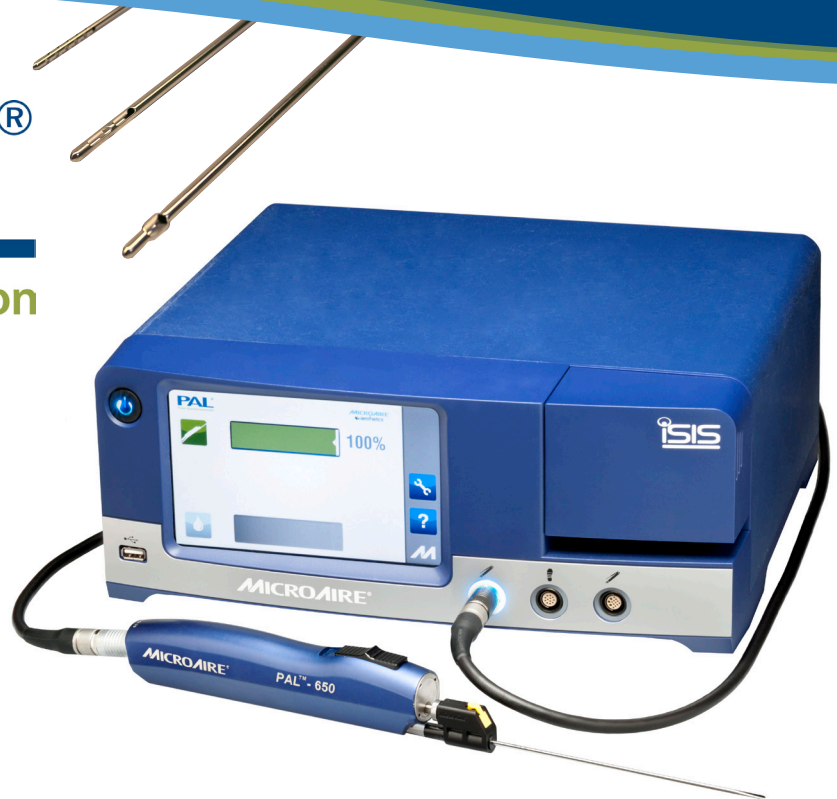

**MICROAIRE®**  
aesthetics  
[microaire.com](http://microaire.com)

**For more information:**

Contact your local MicroAire distributor  
or call MicroAire Customer Service  
+1 434 975 8000.